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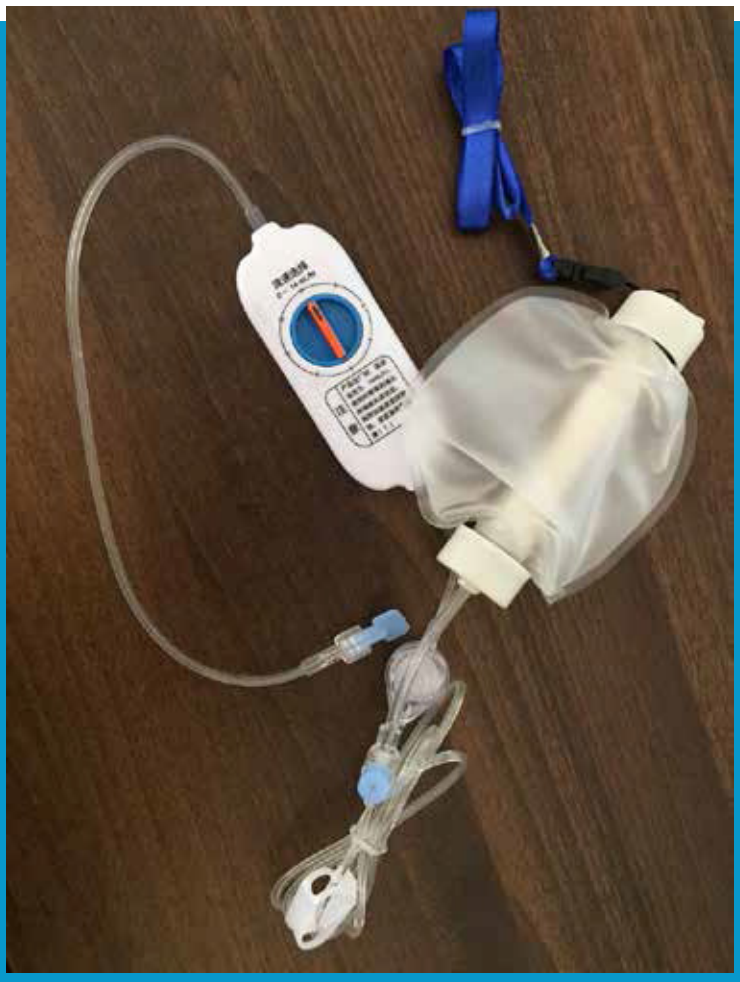
Product's character and usage:

Apply scope:

This product is used for continuous minimal infusion during clinical treatment.

Use instruction:

1 - Take out the infusion pump from the sterile package.

2 - Remove the cap from the add-medicine inlet, clip the stopper closed, infuse the prepared medicine by syringe into the liquid bursa stocker.

3 - If the "filling port" is one-way design, when medicine return is needed, pick out the medicine-return adaptor from the sterilized package, connect the medicine-return adaptor with the syringe, insert into the "add medicine" to return the medicine.

4 - When filling drug, please make sure there is no air in syringe, then infuse the medicine. If there is air in drug balloon, you can turn round and take out all the air (as the drug balloon can exhaust itself. If there is few air, it can exclude after 1~2 hours later.)

5 - After medicine has been added, open the clip, and press the PCA button 1~2 times, for expedite the flow into the pipeline. When the medicine flows from the end, clip the stopper and screw the cover for use next time.

6 - The PCA button is the function button that the patient can use to control the additional medicine when in continuous medicine infusion position. By pressing the PCA button, the patient can add a limited quantity of medicine in respect with doctor's instruction, as he needs.

7 - The flux will be a little faster (within the standard scope) during the 1~2 first hours of use. That is due to the physical characteristic of the silicone bursa.

8 - In the conditions temperature(23± 2) , relative humidity(50± 5)% and atmospheric pressure in 86 KPa~106KPa, the product is made from purified water or distilled water for liquid medium for parallel infusion, mean flow rate for nominal flow rate ± 15%.

9 - Flow rates may vary due to:

1) Fill Volume

When using, please add the liquid according to the nominal volume, too much or too less of fill volume will lead to inaccurate infusion flow.

2) Viscosity and/or drug concentration

The flow rate of the product is calibrated with purified water or distilled water as the medium, infusion of overly viscous liquid will result in a slower flow rate.

3) Temperature

The flow rate of the product is calibrated in the temperature of (23+2) , the flow rate of the product will be fast when the using temperature is on the high side, slow conversely.

4) Atmospheric pressure

The flow rate of the product is calibrated under a standard atmospheric pressure condition, the flow rate of the product will be fast when the using temperature is on the high side, slow conversely.

5) Level infusion

The normal use of the product is level infusion, suspension infusion will lead to a faster flow, slow conversely.

6) Storage

The product should be used timely after the fill volume, the flow rate of the product will be slow after a long time storage.

Construction:

The device is composed of: single valve add-medicine set, drug balloon, medicine liquid filter, bounds flux set, be without), pipeline and conical connectors etc.

Caution:

1 – This device has been sterilized by Ethylene Oxide, it's virgin, no pyrogen.

2 – Check for package integrity before opening. Do not use if package is damaged.

3 – All components are sterile, non-toxic, non-pyrogenic.

4 – For single use only, discard after use.

5 – Please remove the yellow card on the PCA button before use.

6 – Read all instructions before use for patient safety.

7 – Start by infusing the physiological brine that will be needed for dilution of the medicine, then add the medicine, this precaution will avoid un-diluted medicine to flow into the pump line.

8 – Infuse the prescribed dosage of the product, over infusion would effect the result of the cure.

9 – Store in dry place at room temperature, protect from moisture and light.

10- Shelf life is 5 years after manufacturing

、 Contraindication:

1 – This product is not intended for infusion of blood, blood products, lipids, fat emulsions, or Total Parenteral Nutrition (TPN

2 – This product is forbidden for intramuscular injection.